



# Integrity 215

## Steam Sterilizers

Specifically engineered to offer excellent throughput in a space-saving footprint.

Integrity 215 exceeds all standards for capacity—efficient enough for immediate use, productive enough for high-demand environments, and compact enough for most replacement applications. Building on Skytron's history of quality and innovation, the 215 and 215SG are part of a family of sterilizers designed to meet the wide-ranging needs of today's healthcare environments. Whether you're outfitting a new facility, renovating an older one, or replacing outdated technology, Integrity sterilizers offer industry-changing efficiency, capacity, ease of use, and use of space.

### Twice the capacity of similar sized sterilizers

Processes up to six 25-pound trays.

### Easy to use, program, and check cycle status at a glance

Simple touchscreen controls and a large countdown timer and progress bar make daily operation smoother.

### Save on energy, water, and chemicals

Integrated water recirculation reduces water usage by up to 67%\*, without affecting cycle time or increasing footprint.

### Long life and trouble-free operation

The Integrity Series offers fully jacketed stainless steel chambers, chain-driven electric-powered doors, and stainless steel and copper piping standard.

### Superb diagnostics and support

System-wide parts monitoring, on-screen diagnostics, local technicians, and 24/7/365 call center ensure the reliable operation of your equipment.

### Plug-and-play installation

Simplifies replacement of older technology.



# Integrity 215


## Steam Sterilizers

### 24 Preset and custom cycles



#### Processing cycles

- Immediate use 1 (3 minute)
- Immediate use 2 (4 minute)
- Express prevacuum (immediate use 3)
- Prevacuum 270
- Prevacuum 275
- Gravity 1
- Gravity 2
- Liquid cycle (optional)
- Preset cycles validated to AAMI ST-8:2008

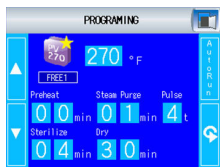


#### Test cycles

- Warm up
- BD
- Cycle
- Leak

#### Custom cycles

12 custom cycles modifiable to your unique needs

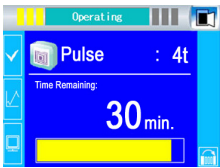

Sterilize quick turnaround devices or even complete orthopedic sets with simple touchscreen control

- Spacious 21½" x 21½" x 38" sterilization chamber
- Processes up to six 25-pound trays (150 pounds per cycle)
- Integrated water recirculation system
- 5.7" LCD full-color touchscreen panel
- Countdown timer and progression bar visible from more than 30'
- 12 standard and 12 fully customizable cycles
- Electric-powered sliding door
- Recessed foot-pedal and touchscreen door control
- Stainless steel and copper piping
- Dual-walled, fully jacketed stainless steel pressure vessel
- Programmed automatic cycle start
- Energy saving daily startup/shutdown timer
- Electronic recording of cycle report data
- Integrated ink-to-paper printer for hard copy cycle reports
- Password protected access to supervisor and service modes
- Compliant to AAMI-ST8:2008 and AAMI-ST79:2010 standards

#### Options

- Integrated steam generator (SG models)
- Load/unload cart and carriage
- Recessed mounting or freestanding cabinet.

### Order Configurations

| Integrity Series        | Door Type | Heating         | Control Power   | Voltage (Opt)        | Mounting             | Equipment #    |                |
|-------------------------|-----------|-----------------|-----------------|----------------------|----------------------|----------------|----------------|
| <b>Integrity™ 215</b>   | Single    | House Steam     | 120V, 60Hz, 10A | NA                   | Recessed             | ASSV-AB09-120R |                |
|                         |           |                 |                 |                      | Freestanding         | ASSV-AB09-120F |                |
| <b>Integrity™ 215SG</b> |           | Steam Generator |                 | 208V, 3pH, 60Hz, 75A | Recessed             | ASSV-AB09-208R |                |
|                         |           |                 |                 |                      | Freestanding         | ASSV-AB09-208F |                |
|                         |           |                 |                 |                      | 480V, 3pH, 60Hz, 33A | Recessed       | ASSV-AB09-480R |
|                         |           |                 |                 |                      |                      | Freestanding   | ASSV-AB09-480F |

\*Water usage is based on comparisons of fully loaded sterilizers with similar-size chambers (20" x 20"), using a 270°F prevacuum cycle with 4-minute sterilize and 20- or 30-minute dry time. Testing showed that the Integrity 215 used substantially less water than the similarly sized STERIS® AMSCO® 400 (35% less) and STERIS CENTURY® (43% less).

Skytron provides safe and efficient healthcare solutions delivered through a local representative network that values long-term relationships and comprehensive service.

Focused on providing the highest quality solutions while delivering the lowest cost of ownership, we prioritize transparency and accountability in our actions.

To learn more, visit our website at [www.skytron.com](http://www.skytron.com)



5085 Corporate Exchange Blvd. SE  
Grand Rapids, MI 49512  
1.800.SKYTRON (759-8766)  
email: [info@skytron.com](mailto:info@skytron.com)  
[www.skytron.com](http://www.skytron.com)



REV4 09-22  
PN L2-010-99