



Flex

Ultrasonic Washer-Disinfecter



Clean Cannulated & Non-Cannulated in the Same Load

The Skytron Flex is an ultrasonic washer-disinfector designed specifically for orthopedic, robotic, cannulated, and non-cannulated surgical instruments.

Cleaner


- Eliminates risks from manual cleaning
- Independent lab tested and validated
- Irrigates with fresh, clean water and detergent
- Supports compliance to AAMI ST79 Standards

Faster

- Triple the capacity of standard ultrasonic washer-disinfectors
- Multi-step cycle cleans and disinfects in less than 40 minutes
- Non-critical water rinse cycles available to shorten cycle time
- Streamlined cleaning process for efficient results

Smarter

- One Skytron Flex is more cost-efficient than multiple ultrasonics and washer-disinfectors
- Simplify operations with sonic only and sonic irrigation cycles
- Innovative technology drives down costs

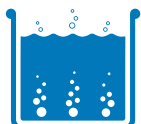

Ultrasonic
Soak



Clean Water
& Detergent
Flush



Ozone
Final
Rinse



125 pound weight capacity includes up to six trays per cycle



Differentiating Features

- Unsurpassed throughput
- Easy to use
- No proprietary trays-use what you currently have
- Fresh water used every cycle
- Cleans orthopedic, robotic, cannulated and non-cannulated instruments simultaneously

Specifications

Dimensions	52" W x 41" D x 42" H
Weight Capacity	125 lbs
HF Output Power	2,000 watts
Frequency	40 kHz with a sweep of ± 2 kHz
Washer Supply	480V/3 phase 50A, 60 Hz
Watts per Gallon	48 Watts/Gallon
Final Rinse	Critical Water Final Rinse

Skytron's Solutions



Sterile Processing

Processing Sinks and Stainless Offerings
Washers, Ultrasonics, and Cart Washers
Prep & Pack Tables & Sterilizers



Clinical

Surgical Tables
Booms & Lights
Full Offering of Stainless Products



Architectural & Environmental

Modular Walls
Custom Window and Door Offering
UVC Disinfection Products



Clinical Business Intelligence

Video Integration
Asset Tracking & Supply Chain Solutions
Surgical Displays

Skytron provides safe and efficient healthcare solutions delivered through a local representative network that values long-term relationships and comprehensive service.

Focused on providing the highest quality solutions while delivering the lowest cost of ownership, we prioritize transparency and accountability in our actions.

To learn more, visit our website at www.skytron.com



5085 Corporate Exchange Blvd. SE
Grand Rapids, MI 49512
1.800.SKYTRON (759-8766)
email: info@skytron.com
www.skytron.com



REV3 09-22
PN L2-010-150